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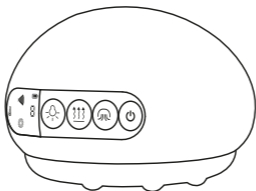
Disclaimer

PRECAUTION

1. If your charging cable or plug is damaged, never disassemble by yourself. Please contact the supplier to request a reparation.
2. Unplug while you finish massaging.
3. Unplug when you are going to clean the product.
4. Never use the parts and accessories that are not recommended.
5. Stop using the product while you feel uncomfortable.
6. Do not pour water or any other liquid onto or into the product.
7. Do not use the product violently.
8. Before using the product, kindly put the filter cotton into the filter. The filter can prevent effectively essential oil from enter into the inner of the product.
9. Wash the filter cotton to reuse.
10. Patients who have varicosity are not recommended to use the product.
11. Do not use the product while you are during the menstruation or pregnancy.
12. Consult a doctor if you have hypoglycemia or hypertension.
13. Never use the product after meal. Kindly take a rest around 30~45 minutes.
14. Take a rest around 30-45 minutes before showering.
15. Scraping ONLY along one direction.
16. To maintain your body healthy, do not scrape every day. Kindly take a rest after scraping.
17. Do not use the product while your body is wet.
18. The product cannot be seen as medicine or medical appliance.
19. People whose skin is sensitive should have a test by smear the essential oil on the wrist or the back of your ear before using the essential oil.

Attention For Environment

1. Please do not use the product in any environment which is extremely hot or humid.
2. The product is prohibited to use while you are driving or walking.
3. It is important to keep in a safe and dry place, away from dust and corrosive odor.

Stop using the product while you feel uncomfortable.

Warning For Scraping


1. Do not scrape on the chest, facial and cervical area.
2. Do not keep scraping the same area.

Attention For After-scraping

1. Scraping helps you boost metabolism. Please drink two glasses of water after scraping, and then take a rest for more than an hour.
3. Scraping will allow your body pores opened, so please be aware of winds after scraping. Taking rest for at least 3 hours before you go to shower. (warm shower would be great.)
4. Take an appropriate rest after scraping. 5-7 days are seen as an standard rest.

Reminder: If you have certain disease, kindly consult a doctor first.

Operation Of Filter Cotton And Filter:

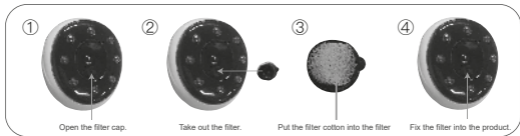The Efficacy and Precautions of Essential Oil

Essential oil efficacy

EN It Soothes emotions, clears the mind, and concentrates; increase the oxygen content of the skin, unblock the pores, and improve the skin.

Instructions for use pure essential oils

Cannot be used directly on the skin. They need to be diluted with base oil or lotion before use. It is recommended to mix ≤ 5 drops of pure essential oils with 10mL base oil.

Precautions

1. Store in a cool place away from light. It is recommended to store below 20 °C/68 °F.
2. Please keep it out of reach of children; it is not recommended for children under 6
3. years old, pregnant women and breastfeeding women.
4. If it gets into eyes accidentally, please wash with plenty of water.
5. Do not try to smell directly near the mouth of the bottle, and do not eat.


For those with particularly sensitive skin for the first time, it is recommended to perform a sensitivity test behind the ear or wrist before use. If there is no allergies, it can be used.

Methods of Scraping

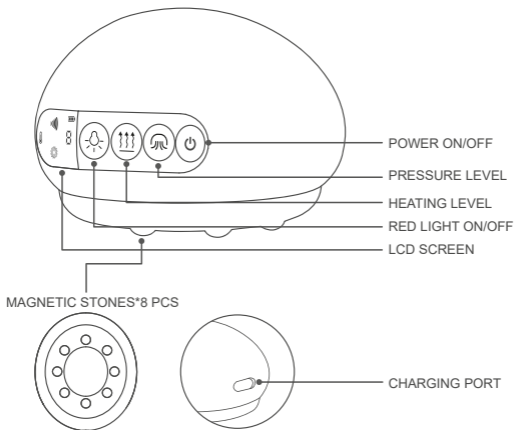1. After applying essential oils, use the following methods.
2. Back: Scrape the back from top to bottom.
3. Sole of feet: Massage from heel to toe.
4. Abdomen: Scrape from top to bottom, or operate in a circular motion.
5. Upper limbs: scrape from top to bottom.
6. Legs: Scrape the legs up from the calf.
7. When using, you should scrape in one direction, not up and down or left and right.

Product Feature

1. The product is hand-held.
2. The outer appearance is ergonomic design.
3. Integrate and multifunctional.
4. Magnetic stones greatly boost your metabolism and blood circulation.

How to Use

1. Long press  to turn on. The red light is on when the product is at standby.
2. Long press  to turn off the product.
3. Short press  to switch the levels (level 0-level 1-level 2-level 3-level 4-level 5-level 0-cycle in turn), and the pressure icon will flash with a long beep accordingly.
Long press the  to turn off the product.
4. Short press  to switch the temperatures (level 0-level 1-level 2-level 3-level 0- cycle in turn), and the heating icon will flash with a long beep accordingly.
Long press  to turn off heating function.
5. Short press  to turn on or off the red light.
 - a. Default Mode: red light is on.
 - b. Pressure Mode: level 1-level 2-level 3-level 4-level 5
 - c. Heating Mode: level 1-level 2-level 3
6. Low Heating: 38°C±3°C/100.4°F±5.4°F
7. Medium Heating: 41°C±3°C/105.8°F±5.4°F
8. High Heating: 44°C±3°C/111.2°F±5.4°F
9. As the product is under charging, it will stop working.
10. Do not use the product while it's under charging.

How to Charge

1. This product comes with a lithium battery. When the power is low, please connect the external power supply accompanied with the product to charge the product;
2. Please use a DC5V adapter, plug the USB end of the charging cable that comes with the product into the adapter, and plug the other end into the USB-C port of the product. The LCD will display the charging status accordingly.
3. To charge the product, you should use a DC5V/1A adapter, or you can also use a computer USB port to charge the product.
4. When the product is charging, the product can't be used.

Battery Removal Instructions

Battery removal step:

1. Cut off the external power supply and remove the product shell;
2. Cut off the two power cords connected to the battery in turn;
3. Remove the battery and use insulator glue to stick to the break part of the power cord;
4. Please put the waste battery in a fixed collection point;
5. Please discard it in accordance with local environmental laws and guidelines.

Attention:

1. The battery must be taken out of the product before discarding the product;
2. When removing the battery, the power must be cut off;
3. The battery should be disposed of safely;
4. Do not disassemble by non-professionals.

Product Parameter

Product Name	Electric chinese bubble
Model No.	AirBubble
Rated voltage	3.7V
Rated Power	4.4w
Standby Current after Shutdown	<20uA
Noise	<65dB
Low Heating	38°C±3°C/100.4°F±5.4°F
Medium Heating	41°C±3°C/105.8°F±5.4°F
High Heating	44°C±3°C/ 111.2°F±5.4°F
Smart Timer	10min
Battery Capacity	2200mAh

Troubleshooting

Follow the instruction below to solve the problems :

EN

Problems	Possible Reasons	Solutions
Failure to turn on the product.	The battery is low.	Charge the product.
	The product is charging.	Use the product after fully charged.
Failure to charge the product.	The plug is not connected well.	Reconnect the plug.

Clearance and maintenance

1. When cleaning and maintaining the product, please turn off the main unit and use wet wipes to clean the product. When it is particularly dirty, you can use a little medical alcohol (75% concentration) for disinfection.
2. For daily maintenance, please wipe it with a soft cloth with a neutral detergent. Do not use corrosive detergents. Do not soak or rinse in water.
3. Prevent water or other liquids from entering the product during cleaning in case of any short circuit and electric shock hazards during use.

Storage

Please place this product out of the reach of children, away from humidity, high temperature, direct sunlight and water splashes. When not in use for a long time, turn off the power and keep it charged once every two months.

Disposal Instruction

The battery used in this product is a special rechargeable battery and cannot be used in other products.

1. Don't throw it into the fire to make it overheated;
2. Do not pierce with nails. Do not disassemble and transform;
3. Do not directly connect the positive and negative electrodes of the battery with metal;
4. Do not store the product with necklaces, clips , etc.
5. Do not charge, use and place in high temperature places such as near the stove or under the scorching sun;

The above situations can cause overheat, fire, break-down, etc.

Warranty

If within one year from the date of purchase, this product fails due to a defect in material of workmanship, Habys will repair or replace the product, or supply necessary components free of charge.

This warranty excludes:

- a) Damage caused by accident, abuse, mishandling, or transport;
- b) Units subjected to unauthorized repair;
- c) Improper use of this product;
- d) Damage exceeding the cost of the product;
- e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises;
- f) Failure to provide the dated proof of purchase.

HABYS Sp. z o.o.

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tel.: +48 13 44 62 788, biuro@habys.pl

www.habys.pl